

Blood Collection Monitor




An ISO 9001:2015 & CE Certified Company

Arnav Scientific Equipment Blood Collection monitor. A new generation portable Blood mixer, ideally suited for blood collection in donation room as well as outdoor camps. Inbuilt key features along with microprocessor based controls ensure smooth & gentle mixing processor volume accuracy, precise monitoring of time & flow rate.



SALIENT FEATURES:

- Precise control on volume collection (Microprocessor based system), best in class accuracy (deviation < 1%)
- Automatic Counter Balancing
- Automatic change over on battery mode in case of power failure (Battery Backup 10 hrs)
- Pause facility to temporarily suspend the donation process
- Gentle & smooth motorized clamping of blood bag tube prevents damage to the live blood cells.
- Broad base tray with non-slip surface accommodates all type of blood bags
- Automatic release of the bag when lifted
- Audio Visual Alarm
 - Flow rate (< 20 ml/min or > 180 ml/min)
 - Completion of collection
 - Battery depletion
- Visual Indication & Displays
 - Set volume (ml)
 - Collection Process status (Run / Stop / Pause)
 - Actual weight (gm)
 - Flow Rate (ml/min)
 - Time of collection (min)
 - Battery status

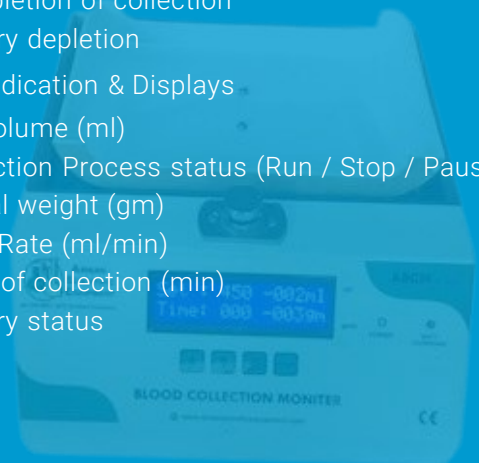

+91 8080007145 / 9973706825

Corp Off.: Gala No. 19, Lusi Industrial Estate,
Behind Hyundai Showroom, Waliv Phata, Sativali Rd.,
Vasai (E), Palghar 401208, Maharashtra, India.

+91 8080007145 / 9973706825 / 7743838033

✉ arnavscientificequipment@gmail.com
sales@arnavscientificequipment.com

🌐 www.arnavscientificequipment.com
www.arnavscientific.com

Blood Collection Monitor



An ISO 9001:2015 & CE Certified Company

Parameter	Unit	ABCM-02
Mixing	Cycle/Minute	12-14
Battery Back Up	Hours	8 to 10
Maximum Volume	ml	Up to 999 net collection
Taring Range	gm	Up to 999
Dimension (W x D x H)	mm	280 x 270 x 255
Net Weight	Kg	5
Electrical	220-240 Volts 50Hz, Single phase	

About Us

Arnav Scientific Equipment is a Manufacturer of the Blood Bank Processing Equipment, Blood Bank Storage Equipment, Blood Bank Collection Equipment. We are manufacturing these products using premium grade raw material that is procured from the authentic vendors of the market. We offer these products at reasonable rates and deliver these within the promised time-frame.

Our Product Range

- Blood Bank Storage Equipment
- Blood Bank Processing Equipment
- Blood Bank Collection Equipment
- Laboratory Equipment
- Blood Bank Equipment
- Blood Bank Refrigerator
- Blood Donor Chair
- Laboratory Refrigerator
- Refrigerated Blood Bank Centrifuge
- Low Temperature Freezer
- Ice Lined Refrigerator
- Stability Chamber

📍 Corp Off.: Gala No. 19, Lusi Industrial Estate,
Behind Hyundai Showroom, Waliv Phata, Sativali Rd.,
Vasai (E), Palghar 401208, Maharashtra, India.

☎ +91 8080007145 / 9973706825 / 7743838033

✉ arnavscientificequipment@gmail.com
sales@arnavscientificequipment.com

🌐 www.arnavscientificequipment.com
www.arnavscientific.com